

# Centri-V 1586

Multi-Purpose Centrifuge

Compare the wide range of Centri-V 1586 features.

The 8-place horizontal, (swinging) rotor allows the processing of multiple samples at once. Gentle and smooth acceleration and deceleration enables you to run fecals with coverslips in place. The unit's innovative design maintains a constant temperature through an innovative vortex ventilation system, maintaining sample integrity.

The Centri-V 1586 allows your hospital to take advantage of horizontal sample separation.

Horizontal separation offers many advantages over traditional fixed-angle separations, including tighter straight-line gel separations to reduce re-spins and re-mixes. The 1586 spins test tubes, up to 16.5 x 128mm, at speeds up to 5,000 RPM, with microprocessor speed control. G-forces up to 4,100g's will produce platelet-poor-plasma in minutes.

## Sample Processing



## Product Highlights

- Process Fecal, Blood, Urine & Microhematocrit samples
- 8-place horizontal (swinging) rotor
- All digital display, touch pad, programmable
- Life-extending, brushless, temperature-controlled motor
- Wide mouth Stainless Steel chamber



# Centri-V 1586

Multi-Purpose Centrifuge

Brushless, temperature-controlled motor extends the operational life of the Centri-V 1586.

A wide mouth Stainless Steel chamber makes loading and clean-up easy and convenient.



## Sample Processing

### Specifications

- 0-99 minute programmable timer
- 800-5000 rpm programmable speed
- G-forces up to 4100g
- Zero-RPM Safety Interlocking metal lid
- Whisper quiet (<70dB)
- Footprint: 18"D x 14"W x 13"H

#### Standard Accessories:

- 8-place horizontal (swinging) head rotor with 15mL tube shields and bumpers

#### Optional Accessories:

- Hematocrit Tube Carrier (to process Microhematocrit samples)
- 4-place 50mL horizontal tube rotor
- 8-place 15mL 45° fixed angle tube rotor
- 24-place 15mL 45° fixed angle tube rotor

Leading Edge provides Veterinary professionals high-performance, competitively priced, reliable Veterinary medical equipment, backed by world-class support and service. Through our Green Initiative, we've made a genuine commitment to support long term environmental sustainability and we're committed to protecting and preserving the environment. As a single source Veterinary medical equipment company, Leading Edge products are available only through a limited and select group of Distributors, who share in our philosophies and standards of excellence, commitment and service.

### Authorized Distributor




#### All digital display with touch pad controls

Programmable speed and time.  
Anti-spill control pad.



#### 8-place horizontal (swinging) rotor

Standard. Includes 15mL tube shields and bumpers.

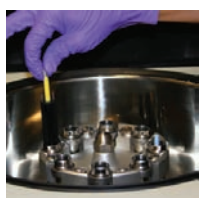

#### Hematocrit Tube Carrier

Allows processing of Microhematocrit samples (optional).