

Axxiom 500

Digital Walk-On Scale

The Axxiom 500 Digital Walk-On Scale combines superior quality with remarkable accuracy and simple operation.


Fast and Accurate.

The Axxiom 500's LCD Digital display provides a rapid weight readout by holding the animal's true weight in less than 10 seconds. Four load cell sensors ensure a 0.2lb weight display accuracy.

Ready for Use in Your Veterinary Hospital.

Designed with the lowest profile available, the Axxiom 500 makes getting animals on and off the scale easier. At only 2" high, this low profile scale also features a durable 22" x 42" stainless steel platform.

Patient Management

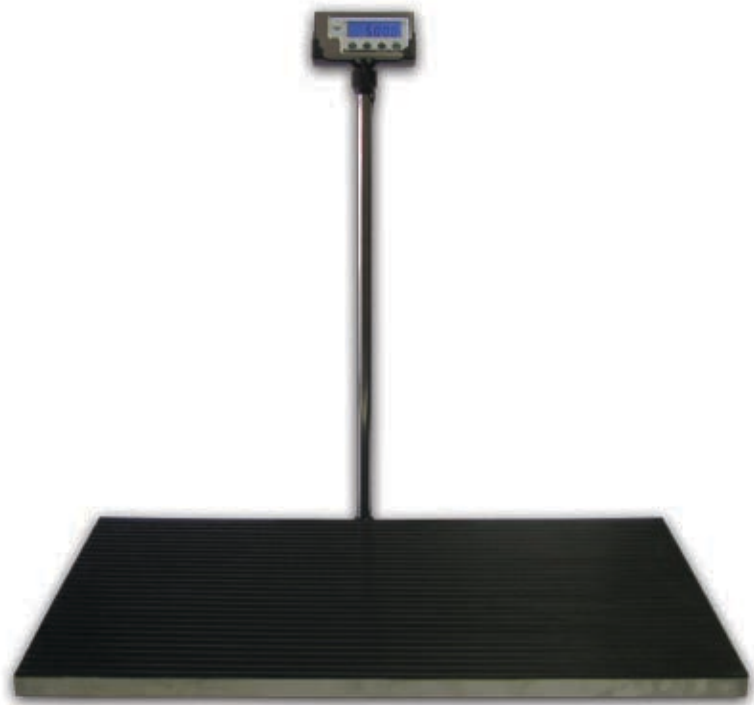

* Floor stand (shown above) optional.



Product Highlights

- Fast, accurate & durable
- Quick weight readings
- Easy to read LCD display
- Precise measurement up to 500 lbs.
- AC or Battery operation
- Auto shut-off saves battery life

Axxiom 500

Digital Walk-On Scale

Patient Management

Simple to Operate.

Four-key push button display including auto zero, tare weight, hold and unit change functions.



Place it anywhere.

The large, easy to read LCD display can be mounted on a wall (with included wall-mount bracket) or on optional Floor or Table-top stands.



Tough, Durable Connections

Solid, tight locking cable connections are durable, resisting moisture and damage. Unlike other scales, they won't come apart when pulled on.



Flexible Configuration

Using either the Table-top or Floor stand accessory allows the digital display to be mounted in the most convenient location.



Level-adjusting Feet

All four feet have rubber pads, allowing use on all floor types and on unlevel surfaces.

Specifications

- 0.1 lb. or 0.1 kg increments
- Accuracy to 0.2lbs
- Stainless Steel Platform 42"L x 22"W x 2H"
- LCD display 6"W x 3H" x 1"D
- Powered by AC adaptor or four (4) "AAA" batteries
- 14' cable
- Four Level-adjustable padded feet

Standard Accessories:

- Display wall mounting bracket
- AC Adapter
- Four (4) "AAA" batteries
- Black textured Rubber Mat

Optional Accessories:

- Table-top stand (19" H)
- Floor stand (36" H)
- Replacement mat

Leading Edge provides Veterinary professionals high-performance, competitively priced, reliable Veterinary medical equipment, backed by world-class support and service. Through our Green Initiative, we've made a genuine commitment to support long term environmental sustainability and we're committed to protecting and preserving the environment. As a single source Veterinary medical equipment company, Leading Edge products are available only through a limited and select group of Distributors, who share in our philosophies and standards of excellence, commitment and service.

Authorized Distributor

