

# Multi-Purpose Table

Lift & Transport Table

## Patient Management

The Leading Edge Multi-Purpose Table is a great value. Big, soft casters deliver a smooth ride that calms animals. Engineered for Heavy Duty use and a long lifespan.

### Truly Multi-Purpose

Veterinarians utilize the Leading Edge Multi-Purpose Table as a treatment table, a vehicle to lift and transport patients in and out of the hospital, perform exams and surgery, and safely lift and transport heavy inventory.

### Reliable, smooth and quiet hydraulic lift

Easy to raise and lower, never needs plugged in and always ready to use. No expensive parts to replace.



## Product Highlights

- Reliable, smooth and quiet hydraulic lift
- Easily lifts and transports 300+ lbs.
- Lift range: 11" - 36"
- Trailing-edge brakes are effortless to engage
- Steering Handle allows safe maneuverability
- Electric lift model also available

# Multi-Purpose Table

Lift & Transport Table

## Patient Management


### Balloon-tread Dual-Locking Casters

Roll smoothly and easily in and out of the clinic. Trailing edge brakes are easy to engage & disengage, won't corrode or damage shoes or pants.



Engineered to provide years of service.

It features a large 22" x 48" stainless steel top with anti-drip edge. It's four corner bumpers protect your walls, doors and door frames.

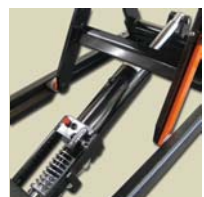

### Maintenance-free Hydraulic Lift

Never needs plugged in. Always ready.



### User-Selectable Casters (Optional)

Provides additional control and safety; makes the Multi-Purpose Table easier and safer to move.



### Qwik-Set Tie Downs (Optional)

Technician's favorite option allows quick and easy patient restraint and release.

## Specifications

- Lift range: 11" - 36"
- Trailing-edge brake tabs are easy to find, effortless to engage and disengage, won't rust or seize and won't tear up shoes or pant legs
- 22" x 48" Stainless Steel top with Anti-drip edge, drain hole and hook for pail
- Four corner bumpers protect your walls, doors and door frames
- Four tie downs, two on each side
- Steering Handle saves your back and provides easy and safe maneuverability without having to lean down to push or pull the table

### Standard Accessories:

- 5" Dual-Locking casters
- Standard tie downs (4)

### Optional Accessories:

- 6" ballon-tread Dual-Locking casters
- User-Selectable Casters
- Qwik-Set Tie Downs
- Restraint Straps
- Electric Lift

Leading Edge provides Veterinary professionals high-performance, competitively priced, reliable Veterinary medical equipment, backed by world-class support and service. Through our Green Initiative, we've made a genuine commitment to support long term environmental sustainability and we're committed to protecting and preserving the environment. As a single source Veterinary medical equipment company, Leading Edge products are available only through a limited and select group of Distributors, who share in our philosophies and standards of excellence, commitment and service.

## Authorized Distributor

