

# Axxiom 2000

Digital Walk-On Scale

## Patient Management


The Axxiom 2000 Digital Walk-On Scale combines superior quality with remarkable accuracy and simple operation.

### Fast and Accurate.

The Axxiom 2000's LCD Digital display provides a rapid weight readout by holding the animal's true weight in less than 10 seconds. Four load cell sensors ensure a 1.0lb weight display accuracy.

### Constructed with Steel to provide years of service.

The Axxiom 2000's large 30" x 59" tread plate surface protects animals and handlers from accidental slips and falls. Featuring easy to walk on ramped edged, and only 2.5" high, this scale accomplishes a multitude of weighing jobs. It fits easily into alleyways for weighing hogs, nursery pigs, sheep, goats or other livestock.

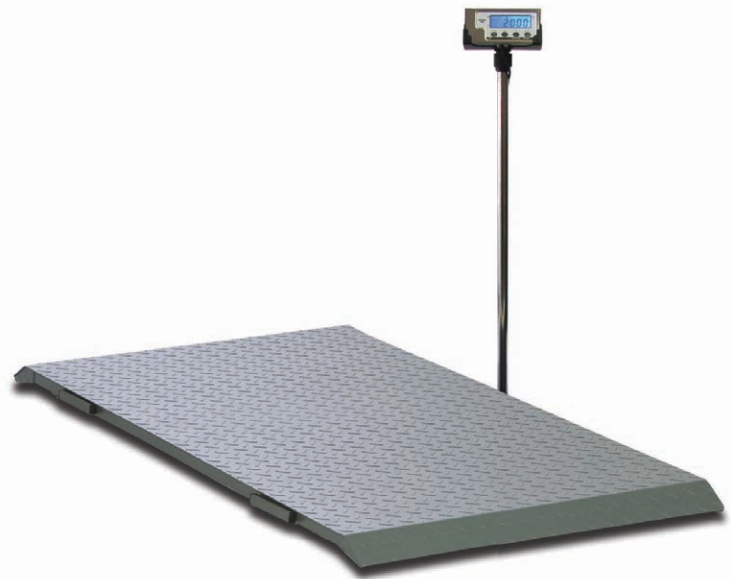

\* Floor stand (shown above) optional.

## Product Highlights

- Fast, accurate & durable
- Quick weight readings
- Easy to read LCD display
- Precise measurement up to 2000 lbs.
- AC or Battery operation
- Auto shut-off saves battery life



# Axxiom 2000

Digital Walk-On Scale

## Patient Management

Simple to Operate.

Four-key push button display including auto zero, tare weight, hold and unit change functions.



Place it anywhere.

The large, easy to read LCD display can be mounted on a wall (with included wall-mount bracket) or on optional Floor or Table-top stands.



### Tough, Durable Connections

Solid, tight locking cable connections are durable, resisting moisture and damage. Unlike other scales, they won't come apart when pulled on.



### Flexible Configuration

Using either the Table-top or Floor stand accessory allows the digital display to be mounted in the most convenient location.



### Collapsible carrying handles

Makes scale easy to relocate.

## Specifications

- 1.0 lb. or 0.5 kg increments
- Accuracy to 1.0lbs
- Steel tread plate platform: 59"L x 30"W x 2.5"H
- LCD display 6"W x 3"H x 1"D
- Powered by AC adaptor or four (4) "AAA" batteries
- 14' cable
- Four Level-adjustable padded feet
- Easy access, collapsible carrying handles makes this scale convenient to relocate
- Operating temperature: 41° - 95°F

### Standard Accessories:

- Display wall-mounting bracket
- AC Adapter
- Four (4) "AAA" batteries

### Optional Accessories:

- Table-top stand (19" H)
- Floor stand (36" H)

Leading Edge provides Veterinary professionals high-performance, competitively priced, reliable Veterinary medical equipment, backed by world-class support and service. Through our Green Initiative, we've made a genuine commitment to support long term environmental sustainability and we're committed to protecting and preserving the environment. As a single source Veterinary medical equipment company, Leading Edge products are available only through a limited and select group of Distributors, who share in our philosophies and standards of excellence, commitment and service.

## Authorized Distributor

