

# Focus-V 1475

## Compound Microscope

### Microscopy

#### Superior Quality at a Competitive Price

The Focus-V 1475's practical design coupled with its built-to-last pathology-grade quality provides more value and features than the competition at a lower price.

#### Perfect for Veterinary Clinicians

The Focus-V 1475 is a compound biological microscope designed for a wide range of uses, ideal for routine veterinary laboratory work and veterinary educational uses.

#### LED Illumination

Powered by an efficient, white and bright LED illumination system, the Focus-V 1475 produces sharp contrast and clear images at any magnification. The LED's cool illumination is perfect for viewing live specimens.



### Product Highlights

- Pathology-grade biological microscope
- Binocular Siedentopf Head - Inclined 30°
- DIN PLAN Objectives: 4X, 10X, 40XR, 100XR (oil)
- Double-layer Mechanical Stage with Side Clamp
- Built-in Variable LED Illumination System

# Focus-V 1475

## Compound Microscope

### Microscopy

#### Specifications

- DIN PLAN objectives: 4X, 10X, 40XR, 100XR(oil)
- 1.25 N.A. ABBE Condenser with Iris Diaphragm
- Binocular Seidentopf Head. Inclined 30°
- Total magnification: 40X -1600X
- Double-layer Mechanical Stage: 5.5" W x 5.5" D
- Hydraulic slide clamp
- Overall dimensions: 14" H x 7" W x 10" D
- Quadruple Reverse Facing Nosepiece
- Variable LED illumination
- Warranty: Limited Lifetime

#### Standard Accessories:

- Wide Field Plane-Scope
- 10x and 16x eyepieces
- Dust Cover
- Filters (Blue, Green, Yellow)
- Immersion Oil


#### Optional Accessories:

- Storage/Transport Case (hard-side)
- Trinocular Head
- Microscope Cleaning Kit
- Immersion Oil

Leading Edge provides Veterinary professionals high-performance, competitively priced, reliable Veterinary medical equipment, backed by world-class support and service. Through our Green Initiative, we've made a genuine commitment to support long term environmental sustainability and we're committed to protecting and preserving the environment. As a single source Veterinary medical equipment company, Leading Edge products are available only through a limited and select group of Distributors, who share in our philosophies and standards of excellence, commitment and service.


### Pathology-Grade Quality Solution

The body of the Focus-V 1475 is equipped with high quality optics to provide for reliable microscopy observation and a broader field of vision. Four flat-field, parcentered and parfocaled DIN PLAN objectives (4X, 10X, 40XR, 100 XR (oil). Compare to Swift in quality, craftsmanship, and performance.

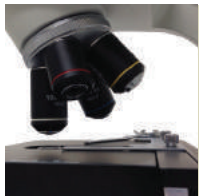

#### Quadruple Reverse Facing Nosepiece

Advanced ball-bearing construction guarantees smooth rotation and precise adjustment.

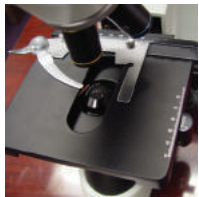

#### Large double layer mechanical stage

Built-in graduated mechanical stage with conveniently positioned coaxial X-Y controls provide smooth-moving lateral and forward/reverse slide positioning.



#### Focus Controls

Coaxial coarse and fine focusing system with durable rack & pinion gearing mechanism for long life. Tension adjustment eliminates stage drift.

Authorized Distributor

ver.070116